

CRUMA1200

Cruma ductless fume hood with the **largest working surface**. Perfect to remove gaseous polluting agents and/or solid particles/aerosols from the working space in a simple, efficient and cost effective way, protecting both the user and the environment. CRUMA 1200 ductless fume hood uses the CRUMA filtration system, without any exterior connection.


PROTECTION OF AND RESPECT FOR THE ENVIRONMENT –The toxic chemical products are not released outside but retained within the filter..

COST SAVINGS AND FAST SET UP –No building work is needed to install ducts to channel the gas outside, which means fewer problems for a laboratory working at full capacity.


FLEXIBLE –It can be used in areas where it is difficult to remove/extract contaminated air, such as from the lower levels in buildings that have a number of floors.

ENERGY SAVINGS –The air that is sucked in is not expelled but recirculated back into the lab after the contaminating substances have been removed. This means that it is not necessary to increase the use of the air-conditioning or heating systems to compensate for the air removed.

CERTIFIED –Manufactured in Spain and certified by an external laboratory according to international standards and applying the criteria of ISO 9001.

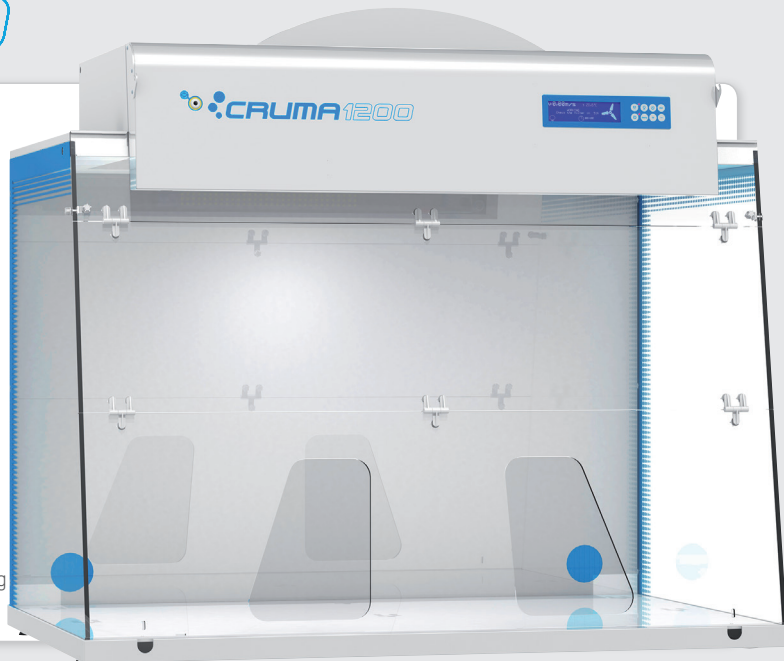

NEW FEATURES

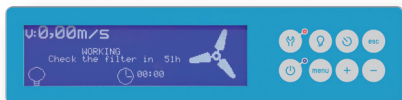

More information on the new LCD display


- ✓ New size 127x34mm display
- ✓ Air speed continuously monitored
- ✓ Type of filter installed, working hours, expiration date and next revision date
- ✓ Open door warning through electronic photocell
- ✓ Countdown timer
- ✓ Clock and calendar

New features and components

- ✓ Initial air flow cycle adequacy and final purge cycle
- ✓ Fault LED
- ✓ Control of air flow through Microprocessor
- ✓ Activated carbon filters with electronic chip
- ✓ LED illumination

New alarms and scheduled warnings

- ✓ Open door warning
- ✓ Open door in off mode warning
- ✓ 60h of filter use warning
- ✓ Next validation warning
- ✓ Few hours of filter life warning
- ✓ Countdown timer warning
- ✓ Expired filter alarm (by hours)
- ✓ Expired filter alarm (by date)
- ✓ Temperature alarm
- ✓ Equipment without filter alarm
- ✓ Low barrier alarm

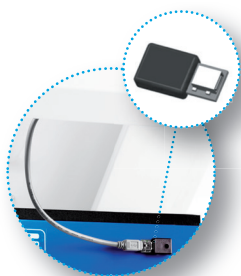

USES

General chemistry involving small volumes of reagents or chemical compounds at ambient/moderate temperature in all types of laboratories:

- ✓ Research laboratories
 - ✓ Quality control laboratories
 - ✓ Clinical and hospital laboratories
 - ✓ University and school laboratories
- ...In general, in any kind of laboratory.

TECHNICAL FEATURES

Number of filtration columns	1
Number of filters	1 to 2
Number of IP44 fans	1
Average volume of treated air	175 m ³ /h
Average face velocity	0,50 m/s
Internal volume of the cabinet	0,692 m ³
Renewals inside the cabinet / min	4,2
Total electrical power consumption	123 W
Voltage-Frequency	220 V - 50 Hz
LED light intensity	36 W / 800 Lux
Noise level	48 dB
Packaging: 100% recycled wooden box with international phytosanitary certificate	Volume 0,64 m ³ Weight 118 Kg

SIZES (mm)

External			Internal		
Width	Depth	Height	Width	Depth	Height
1200	800	1100	1175	760	775

It is not a typographical error,
7 year warranty

Because we are convinced of the quality of our products.



*More information
www.cruma.es

Well done, well shipped.
Our responsible packaging

Wood box 100% recyclable with international phytosanitary certificate.




Do you need help or
technical assistance?

Contact your distributor or call us if you have any questions or need technical support, spare parts, maintenance service...

+34 93 370 61 62

CRUMA1200

SIZES (mm)

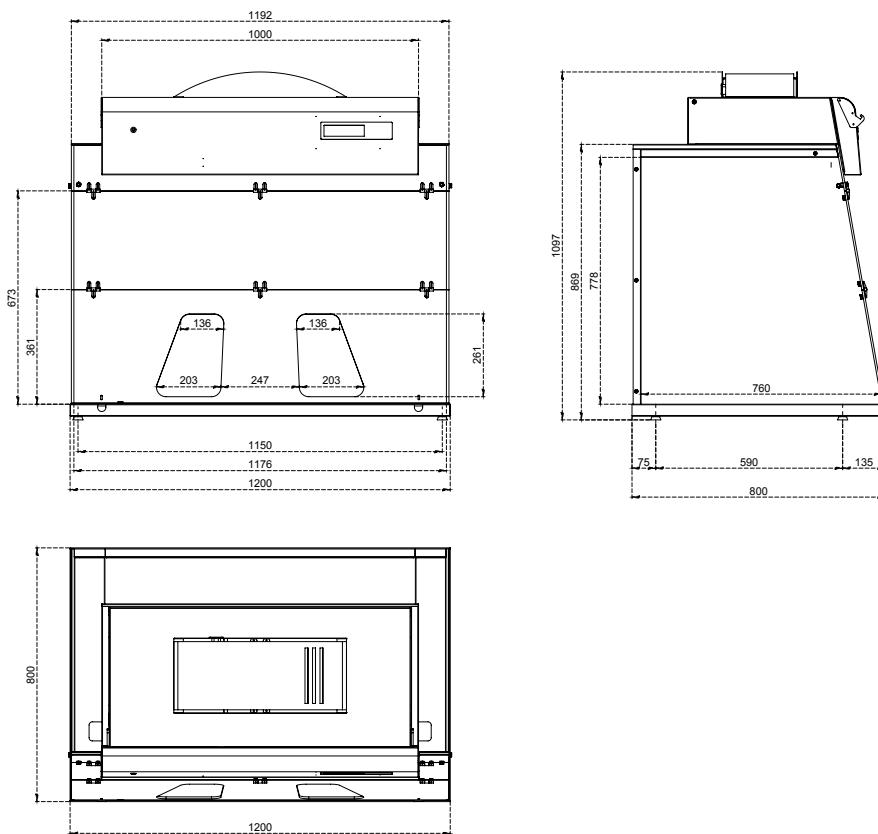

SERIAL EQUIPMENT

Electronic circuit with large format LCD screen	Security levels: level 1 for users and level 2 for maintenance users
Electronic anemometer device	Electronic sensor monitoring continuously air face velocity
Photocell sensor device for open door detection	Electric device with open door alarm
Electronic control device for filters replacement	The filters incorporate a microchip with miniUSB connection that identifies the type of filter installed, the expiry date and the serial no.
Illumination	96 LED Tube high light intensity and low power consumption - 16 Watts / 700 Lux
Sampling system to analyze the filtered air at the exhaust	To detect the level of filter saturation
60 hours alarm	Countdown timer according to French NF X 15-211:2009
Electronic cronometre with audible alarm	To program the work inside the fume hood
Clock and calendar	Display of date and time
Working surface 1	Spill retention tray (2-10 liters) with working surface in white tempered glass
64 Prefilter	64 class pre-filtering blanket of synthetic biofibres (according to EN-779) for the retention of atmospheric dust
Cable entry holes (2)	Access to the rear wall for cables and / or services entry
Chemical Listing	Guide of retained products by type of filter
Warranty	7 years

OPTIONAL EQUIPMENT

Movilair	Stand with wheels and internal tray in Epoxy coated steel
Tubular steel stand	Support stand in Epoxy coated steel
Working surface 2	Spill retention tray (2-10 liters) with working surface in phenolic resin
Working surface 3	Spill retention tray (2-10 liters) with working surface in inox steel
Transparent rear back pannel	Transparent polymethylmethacrylate rear pannel 8 mm thick (light transmission of 93%). Ideal for teaching sessions
Voltage / Frequency	125 V / 50 Hz
Filter test kit	Dräger pump with reactive colorimetric tubes (pack 10u)

MAIN STRUCTURE

Metal parts: base frame, rear wall and head	1.2 mm galvanized coated steel with anti acid polymer resin powder heat-hardened at 200 °C
Front and side panels	Transparent polymethylmethacrylate 6 mm thick (light transmission of 93%)



CRUMA 1200

FILTER TYPES

Type A	For organic vapors such as ketones, ethers, alcohols, xylenes ... Eventually it can be used for inorganic acids, but only if used in small quantities because this activated carbon is not impregnated and the excess of acid vapors could saturate it quickly.	Type K	For NH3 vapors and amines ; also good for other organic compounds. Carbon with metal salt complexes impregnation.
Type BE	For inorganic acid vapors as H2SO4, HCl, HNO3, and volatile sulfur compounds such as H2S, SO3, ... It can be used with organic vapors because the activated carbon incorporates impregnation of metal compounds and neutralizing salts. It is also suitable to filter organic and inorganic compounds when they are in similar proportions.	Type ABEK	Mixed type to be used when the ratios between organic, inorganic and NH3/amines are similar.
Type F	For formaldehyde vapors and derivatives; also good for other organic compounds. Carbon impregnated with Cu leads, so that it should never be used with inorganic acid vapors.	Type D	HEPA H-14 filter (High Efficiency Particulate Air, according to EN-1822: 1998) for filtering dust and smoke particles.

MODULAR FILTRATION COLUMN FOR GASES AND PARTICLES (according to NFX 15-211:2009)

CLASS 2

Type G Handling of liquid compounds/products		Type GS Handling of liquid and particles compounds/products	
--	---	---	---

FILTRATION COLUMN FOR POWDERS

Type D
Handling of powder compounds



Fan




Molecular Filter
HEPA-H14 Filter

ACCORDING TO STANDARDS

Cabinets / Fume Hoods	AFNOR NF X 15-211:2009 (France) BS 7258: 94 (UK) BS 7989: 2001 (UK) ANSI/ASHRAE 110-1995 (USA)
Filters	CSA Z 316.5-94 (Canada) EN-779: 1996 (HEPA & ULPA Filters) EN-1822:1998 (HEPA & ULPA Filters) EN-141:2001 (Gas Filters)
Quality	UNE EN ISO 9001:2008

Filters with microchip and USB connection,
intelligent filtering

Information storage and interaction
with the alarm system.

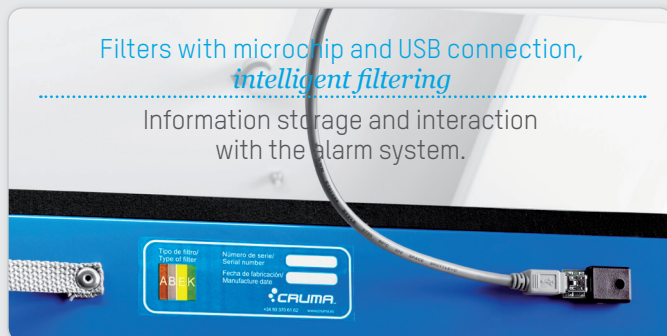


we recognise our responsibility and dependence towards a healthy environment and, therefore, donate 1% of our annual sales to environmental organisations around the world.

A. Canals

our **1%** engagement

CRUMA
PLUG&PLAY RANGE

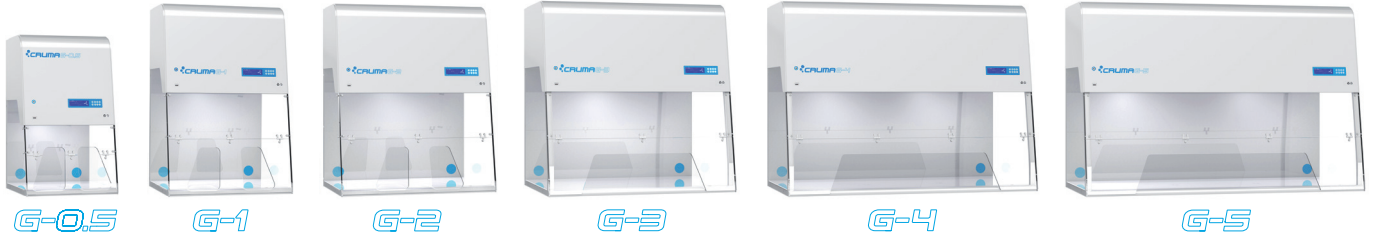

CRUMA
CLASSIC RANGE



CRUMA
POWDER CABINETS




CRUMA
PCR CABINET



CRUMA
LAMINAR FLOW CABINETS



CRUMA
BIOSAFETY CABINETS

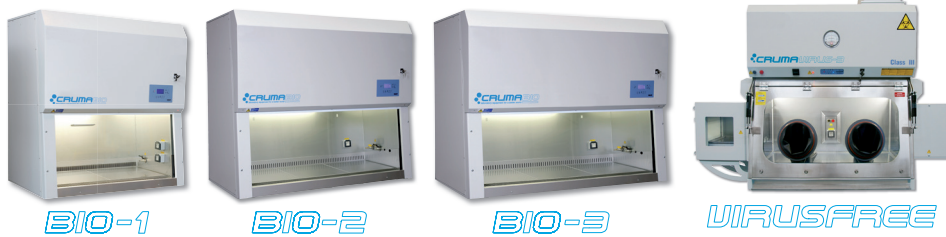

CRUMA
CO₂ INCUBATOR



DISTRIBUTED BY: